

## Instrument controller software run summary:

**Filename and data path:** C:\Agilent Technologies\Data\2022 09 29\12-24-47\2022 09 29 12H 24M.raw

**Created:** Thursday, September 29, 2022 12:40:21 PM

**Number of capillaries:** 2

**Array serial number:** 030918-02USFS

**Effect length:** 22 cm

**Array usage count:** 1577

**Instrument type:** 5300 Fragment Analyzer

**Instrument controller software version:** 3.1.0.12

**Device serial number:** 2551

## Method Information

**Method name:** DNF-474-22 - HS NGS Fragment 1-6000bp.mthds

**Gel prime:** No

**Full conditioning:** Yes

**Gel prime to buffer:** No

**Gel selection:** Gel 1

**Perform prerun:** 7.0 kV, 30 sec.

**Rinse:** No

**Marker 1:** No

**Rinse:** Tray: Marker, Row: A, Dip count: 1

**Sample injection:** 4.5 kV, 30 sec.


**Separation:** 7.0 kV, 25.0 min.

**Tray name:** Tray-1



**Analysis mode:** NGS


## Notes

## Gel Image

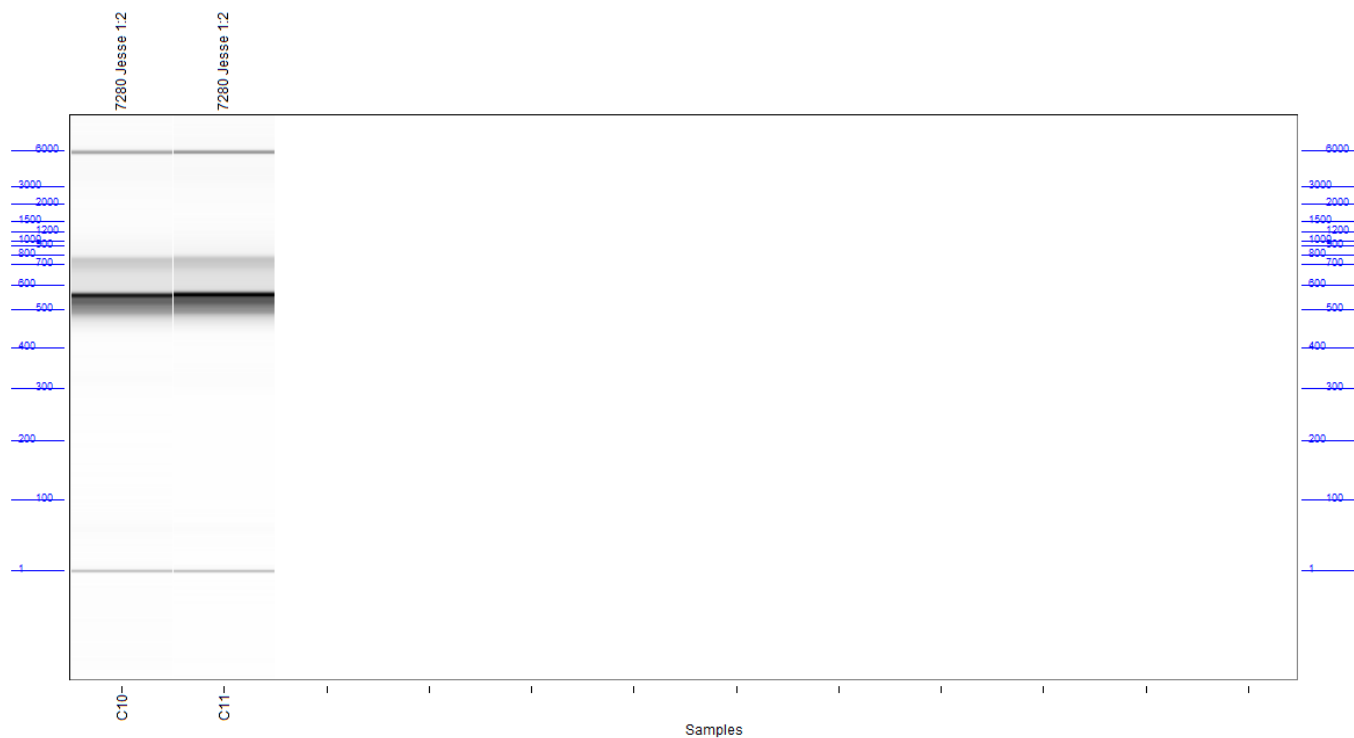

**Filename and data path:** C:\Agilent Technologies\Data\2022 09 29\12-24-47\2022 09 29 12H 24M.raw

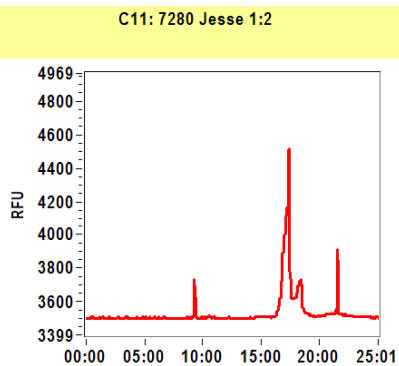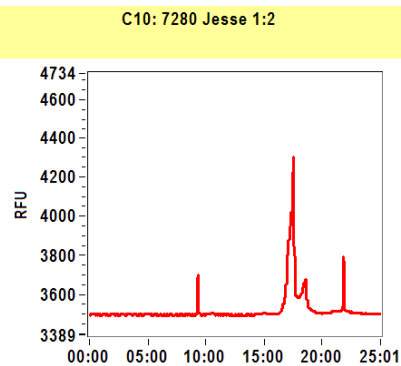


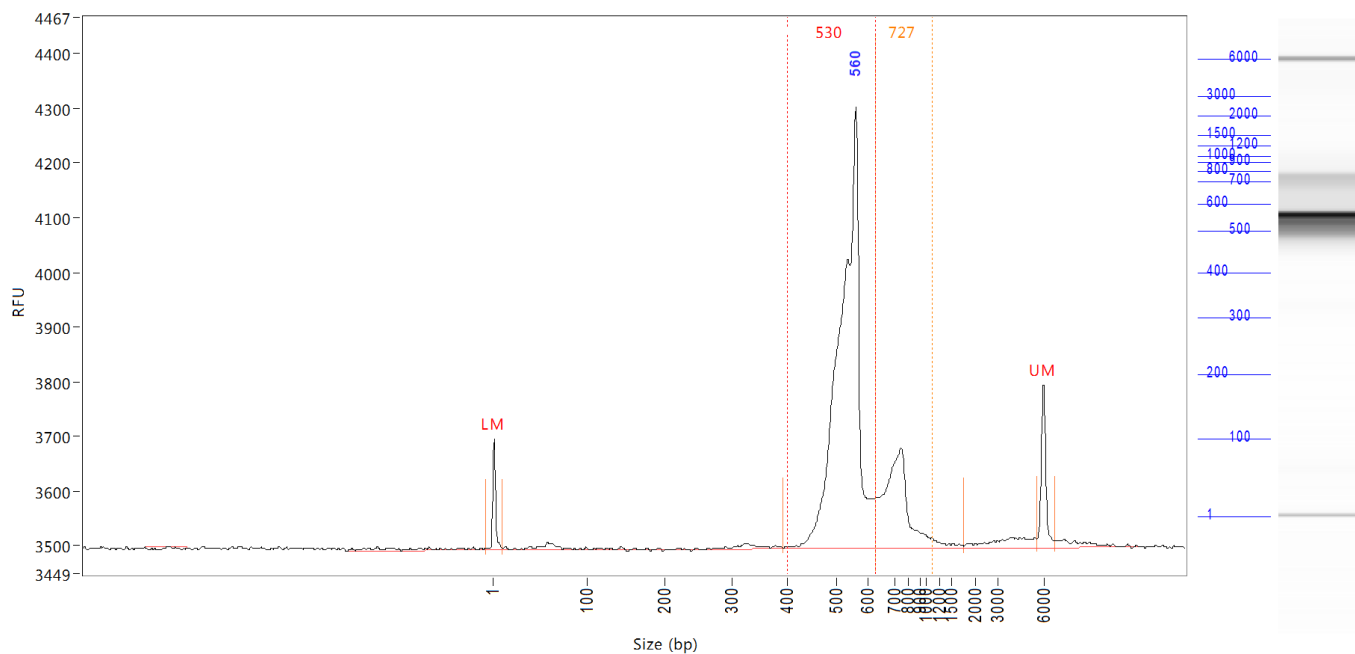**Sample:** 7280 Jesse 1:2**Well location:** C10**Created:** Thursday, September 29, 2022 12:40:21 PM

Peak	Size	Concentration	From	To	Average size	CV%	RFU	Corrected peak area
	(bp)	(ng/uL)	(bp)	(bp)	(bp)			
1	1 (LM)	0.0062	0	10	1	365.85	204	1.677
2	560	1.3054	393	1766	578	21.39	809	29.589
3	6000 (UM)	0.0062	5578	6696	5965	3.13	297	1.688
	TIC:	1.3054	ng/uL					
	TIM:	3.7167	nmole/L					
	Total	1.4097	ng/uL					
	concentration:							

Smear Analysis	400 bp to 625 bp	1.0259 ng/ul	72.8 %Total	3.1833 nmole/L	530 Avg. Size (bp)	6.95 %CV
	625 bp to 1100 bp	0.2704 ng/ul	19.2 %Total	0.6123 nmole/L	727 Avg. Size (bp)	12.07 %CV

Sample peak width (sec): 50    Sample min peak height: 25    Sample baseline V to V?: Y    Sample baseline V to V points: 3  
 Sample filter: Binomial    Number of points for filter: 3    Sample start region (min): 0    Sample end region (min): 25  
 Manual baseline start (min): 6    Manual baseline end (min): 24  
 Marker peak width (sec): 3    Marker min peak height: 200    Marker baseline V to V?: Y    Marker baseline V to V points: 3  
 Lower marker selection: First peak > 200 RFU    Upper marker selection: Last peak > 200 RFU  
 Ladder size (bp): 1, 100, 200, 300, 400, 500, 600, 700, 800, 900, 1000, 1200, 1500, 2000, 3000, 6000  
 Quantification using: Ladder    Final concentration (ng/uL): 0.0830    Dilution factor: 12.0

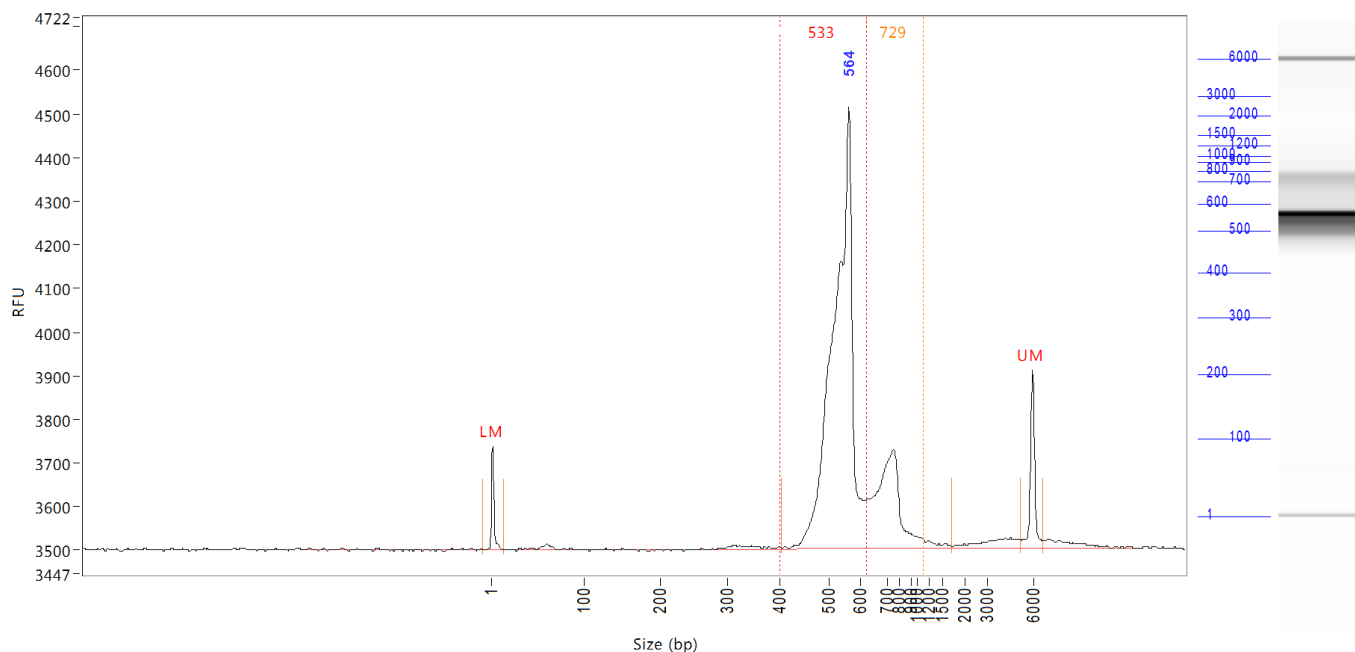**Sample:** 7280 Jesse 1:2**Well location:** C11**Created:** Thursday, September 29, 2022 12:40:21 PM

Peak	Size	Concentration	From	To	Average size	CV%	RFU	Corrected peak area
	(bp)	(ng/uL)	(bp)	(bp)	(bp)			
1	1 (LM)	0.0062	0	13	0	414.43	237	1.916
2	564	1.4067	402	1719	581	21.18	1014	36.417
3	6000 (UM)	0.0072	5203	6696	5946	3.88	407	2.229
	TIC:	1.4067	ng/uL					
	TIM:	3.9843	nmole/L					
	Total	1.4845	ng/uL					
	concentration:							

Smear Analysis	400 bp to 625 bp	1.1021 ng/ul	74.2 %Total	3.4017 nmole/L	533 Avg. Size (bp)	6.81 %CV
	625 bp to 1100 bp	0.2943 ng/ul	19.8 %Total	0.6647 nmole/L	729 Avg. Size (bp)	11.98 %CV

Sample peak width (sec): 50    Sample min peak height: 25    Sample baseline V to V?: Y    Sample baseline V to V points: 3  
 Sample filter: Binomial    Number of points for filter: 3    Sample start region (min): 0    Sample end region (min): 25  
 Manual baseline start (min): 5    Manual baseline end (min): 24  
 Marker peak width (sec): 3    Marker min peak height: 200    Marker baseline V to V?: Y    Marker baseline V to V points: 3  
 Lower marker selection: First peak > 200 RFU    Upper marker selection: Last peak > 200 RFU  
 Ladder size (bp): 1, 100, 200, 300, 400, 500, 600, 700, 800, 900, 1000, 1200, 1500, 2000, 3000, 6000  
 Quantification using: Ladder    Final concentration (ng/uL): 0.0830    Dilution factor: 12.0

**Sample:** NGS Ladder**Well location:** C12**Created:** Thursday, September 29, 2022 12:40:21 PM**Fit type:** Point to point

Calibration curve

